


A1
FACEPANEL

Face Panel
Cut 1
Clear 20mil Vinyl


B1
FACEPANEL

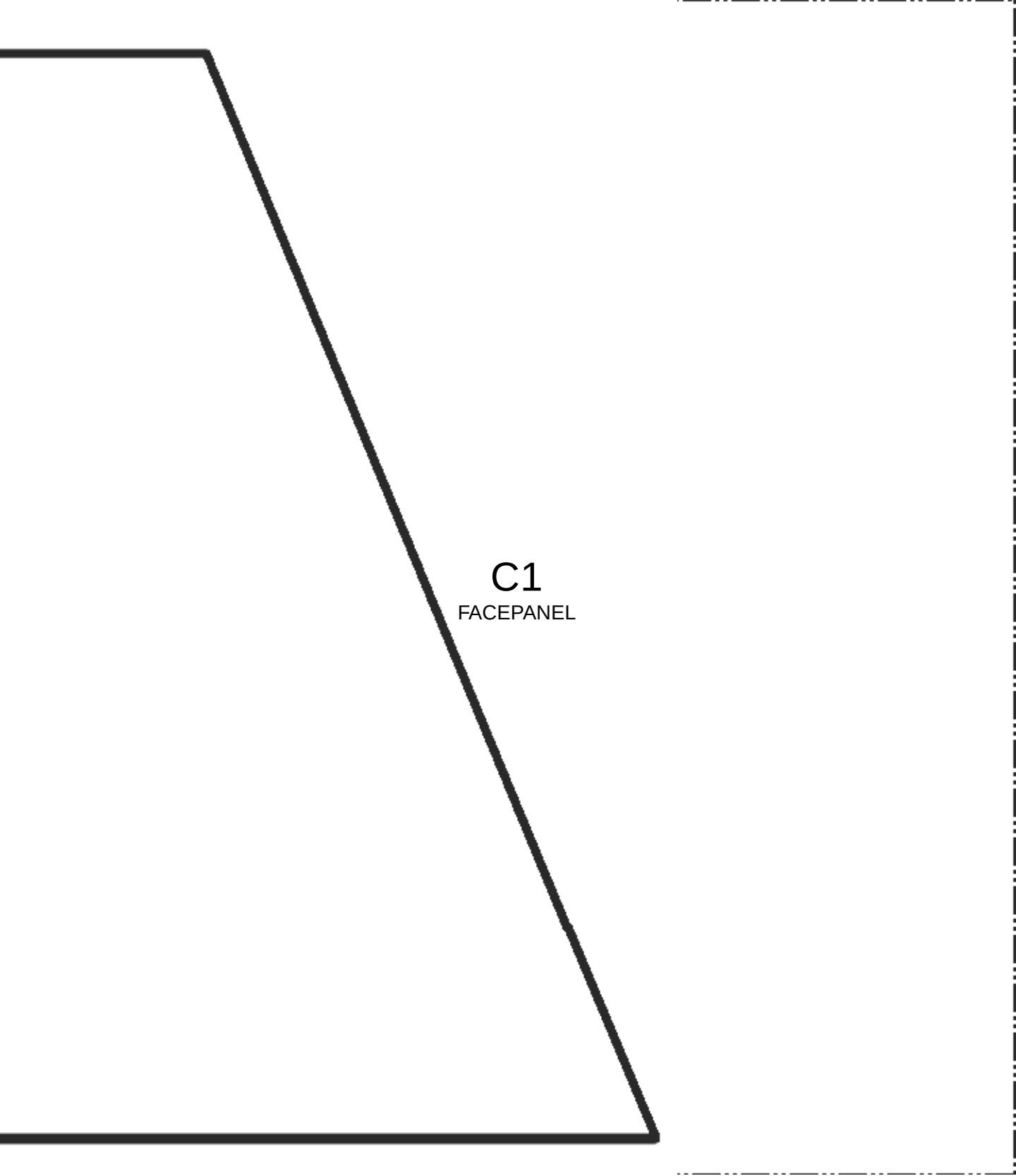

C1
FACEPANEL


1"

A1
Hood Skirt

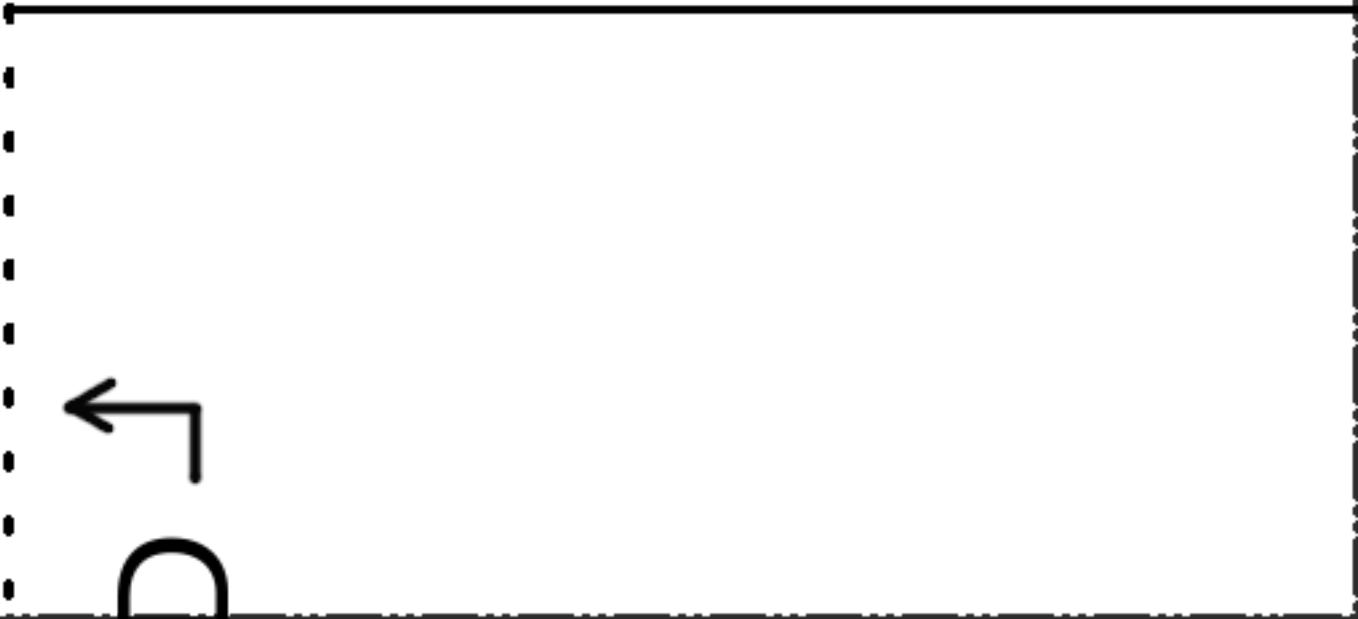

↑ PLACE ON FOLIO

A2
Hood Skirt

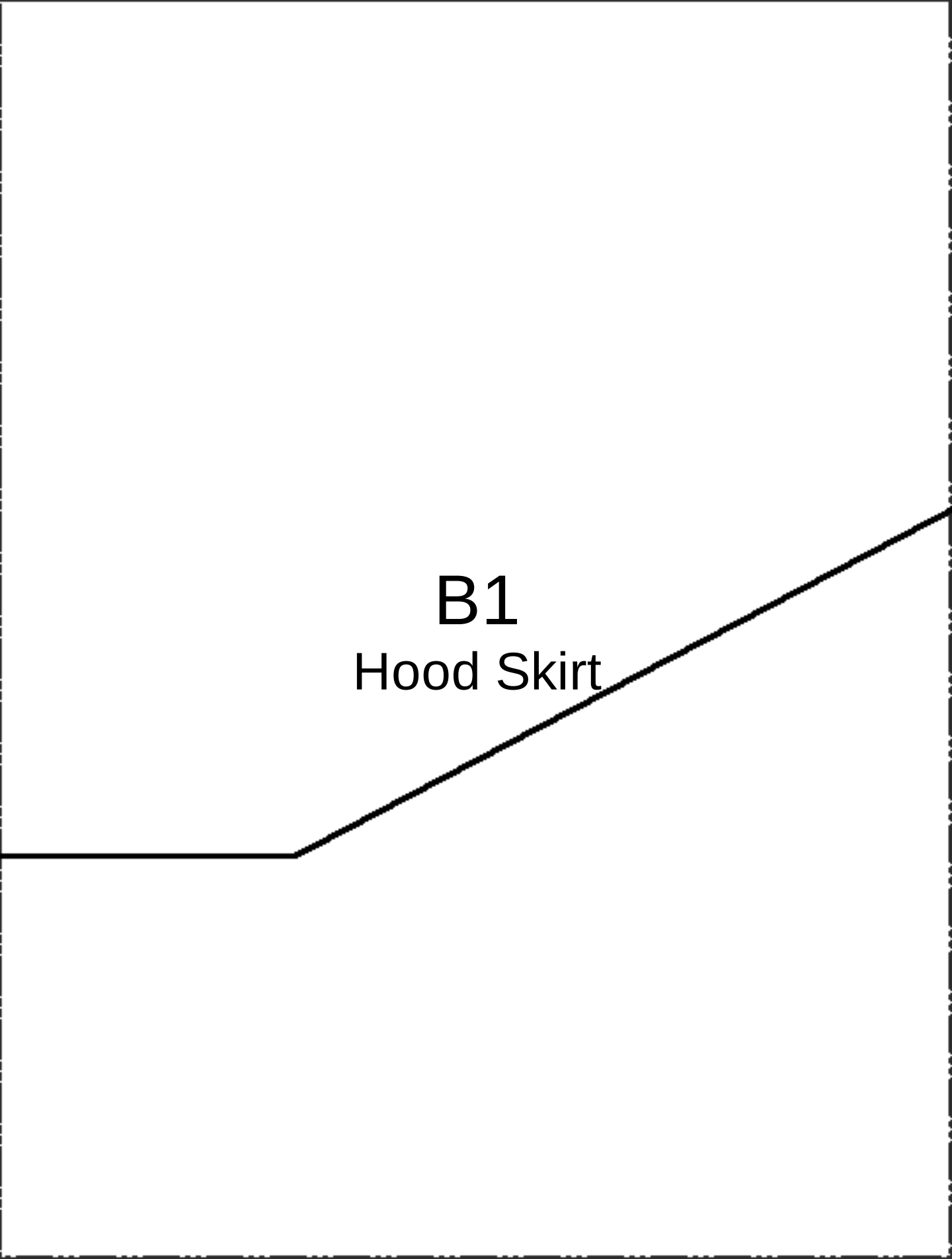



The image shows a line graph on a white background. The graph consists of a solid black line that starts at the left edge of the frame, moves horizontally to the right, and then turns upward at an angle towards the right edge. The text 'B1 Hood Skirt' is centered in the area between the horizontal and the upward-sloping segments of the line.

B1
Hood Skirt

1"

B2
Hood Skirt

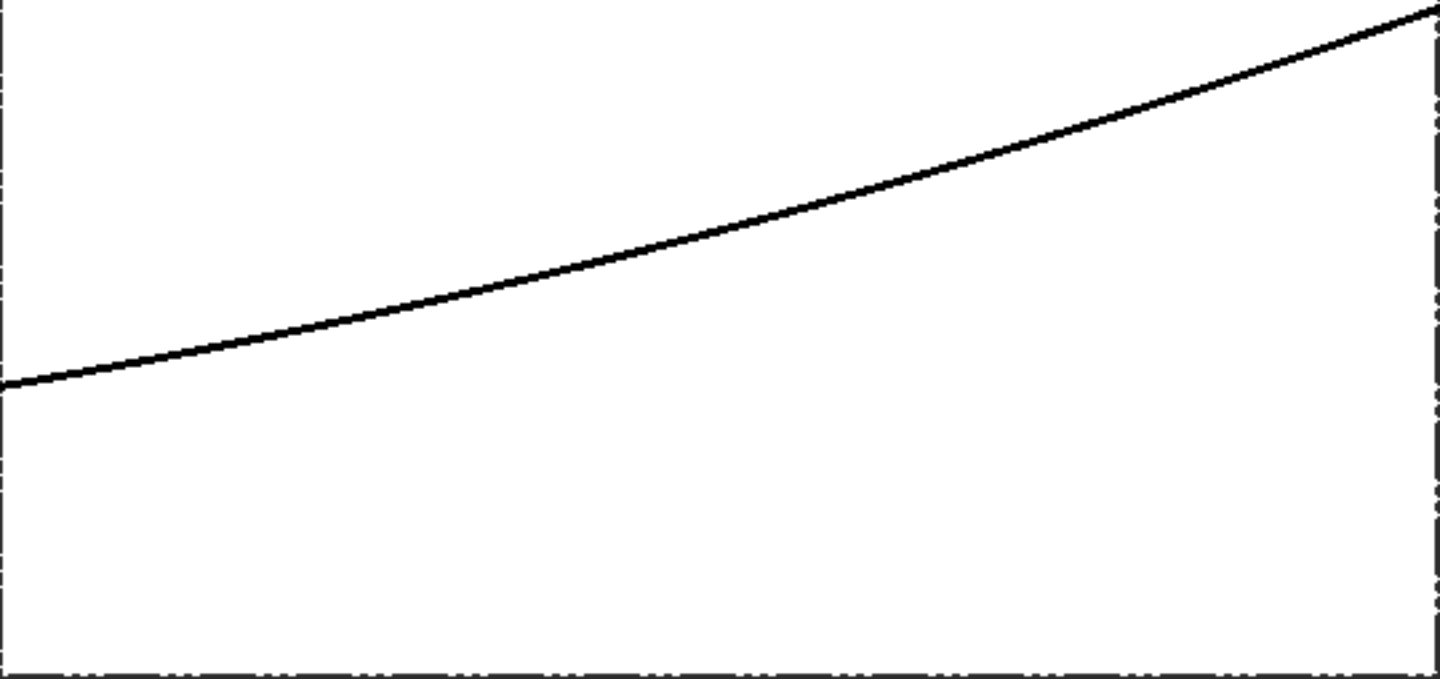

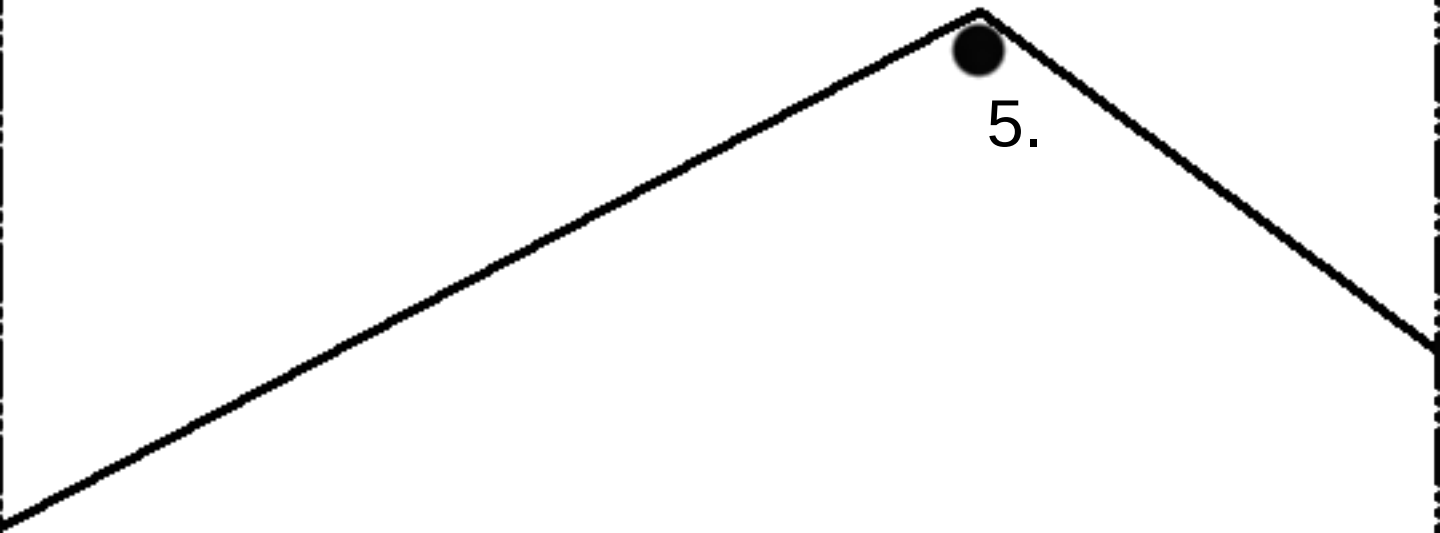


C1
Hood Skirt

HOOD SKIRT
Cut 1 on fold
4mil plastic sheeting




C2

Hood Skirt

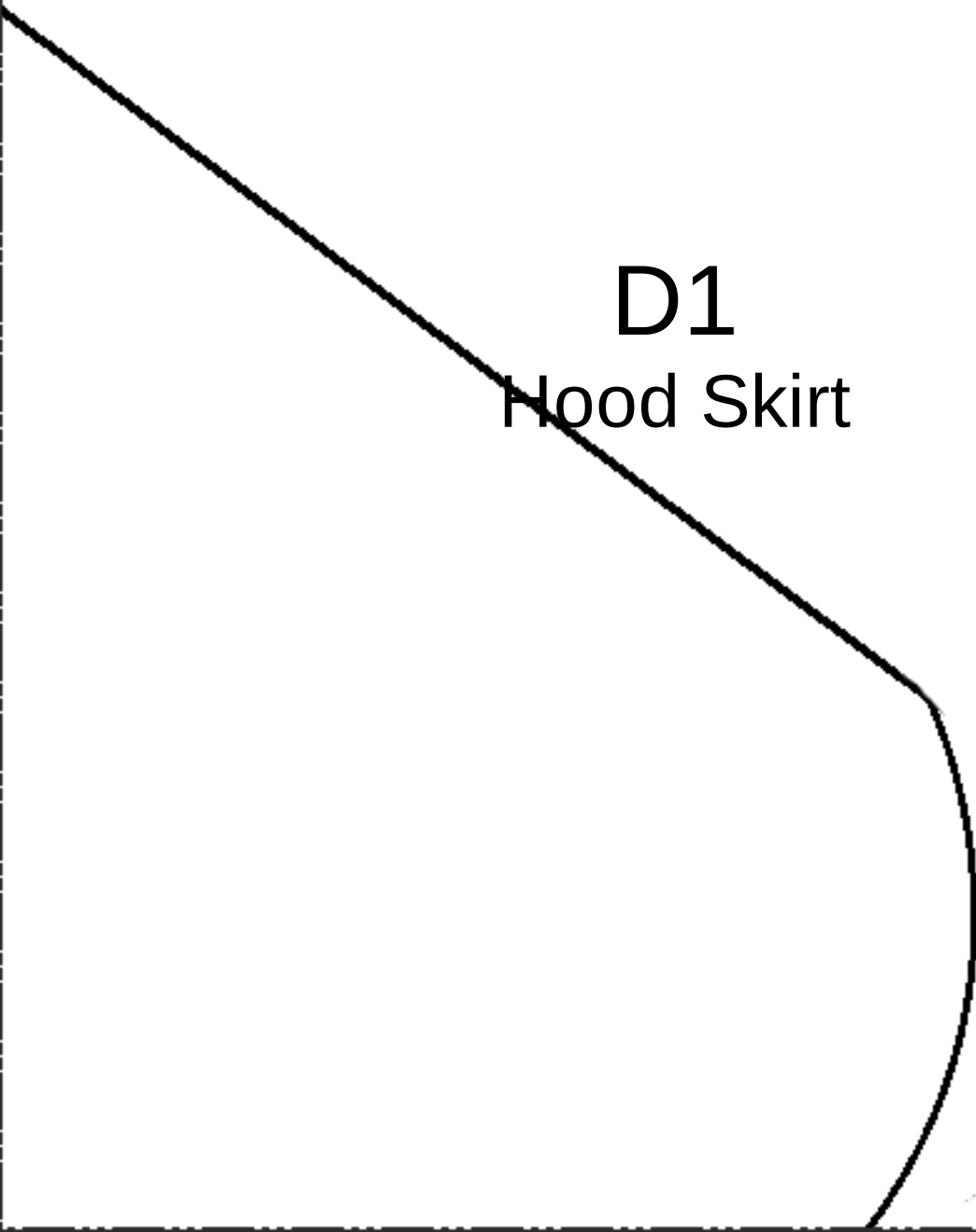

D1

Hood Skirt

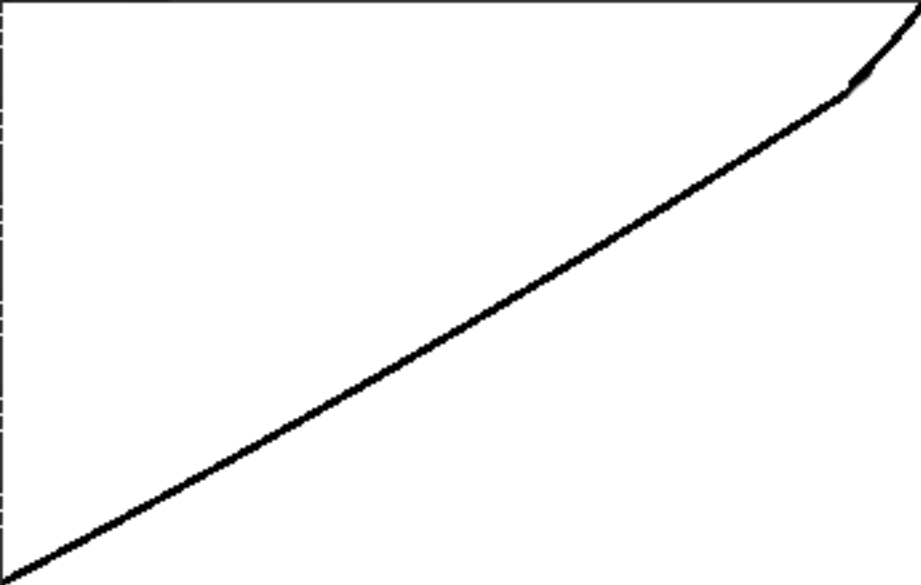


D2
Hood Skirt

1"

3.

A1
Hood Body

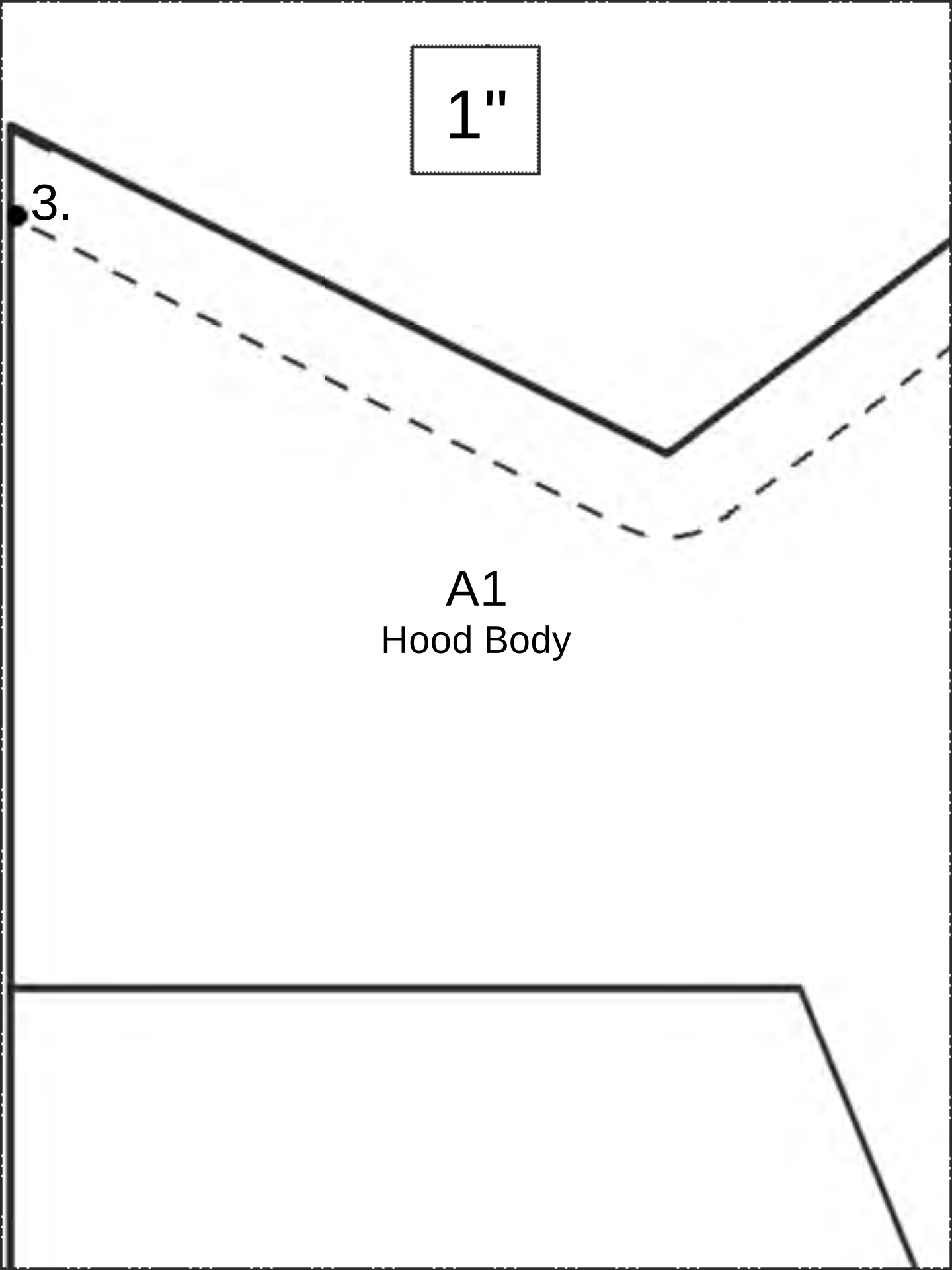


A2
Hood Body


FOLD ↖

↑ PLACE ON I

A3
Hood Body

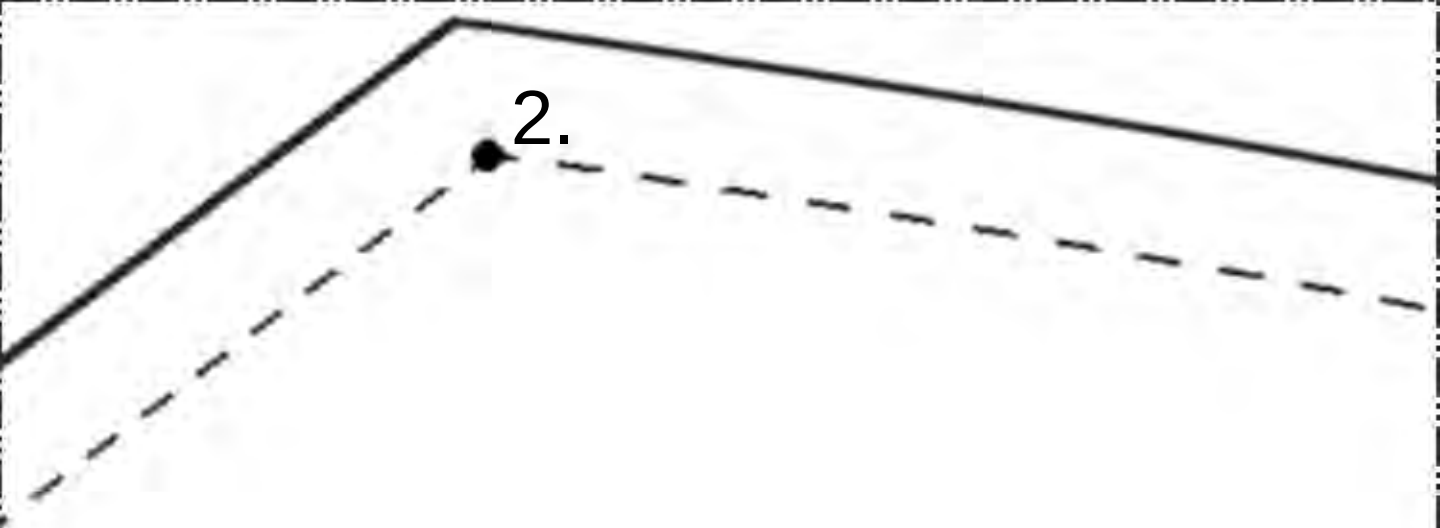

B1
Hood Body

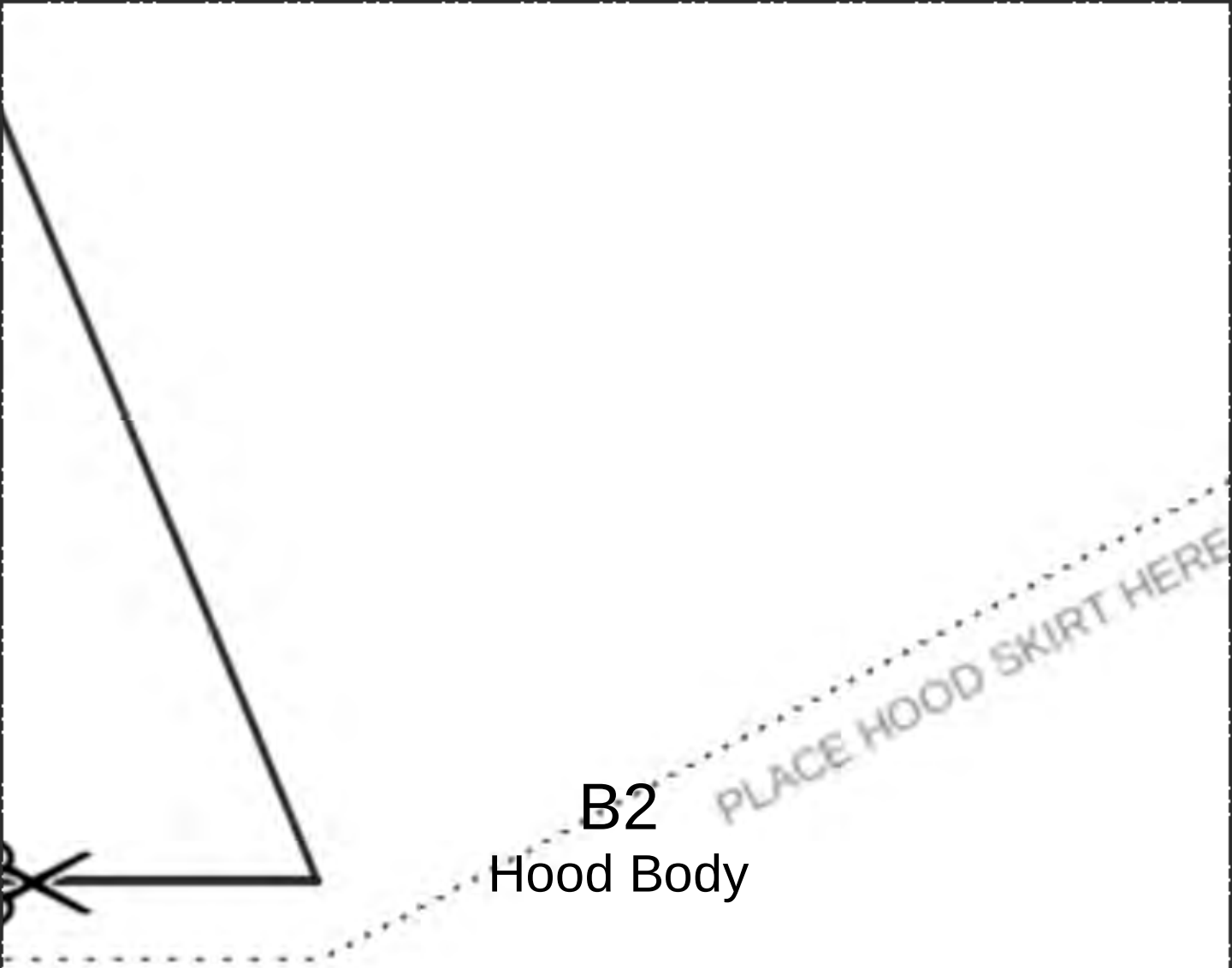

B2

Hood Body


PLACE HOOD SKIRT HERE

HOOD BODY


Cut 1 on fold
6mil plastic sheeting

B3
Hood Body



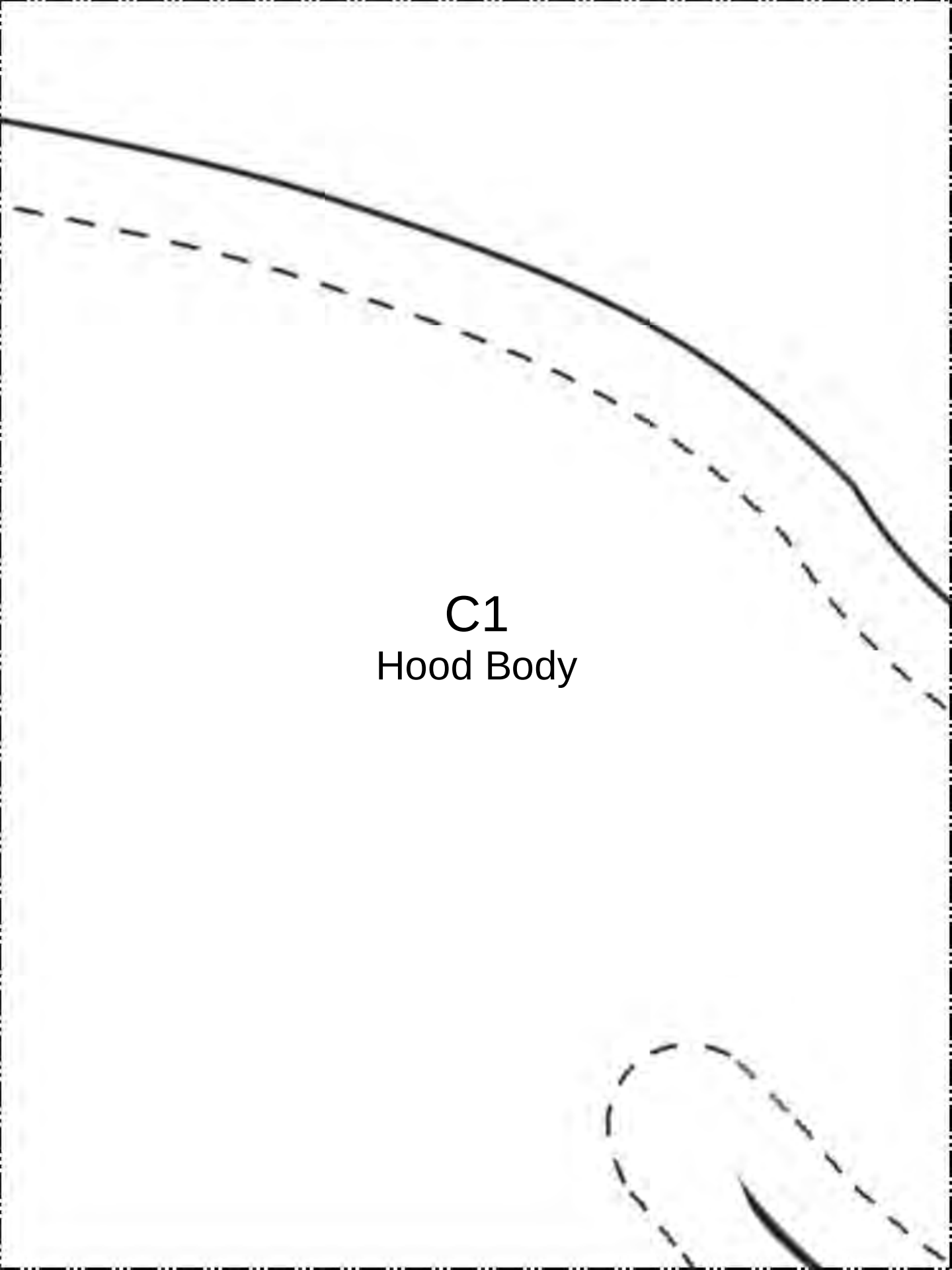


C1
Hood Body

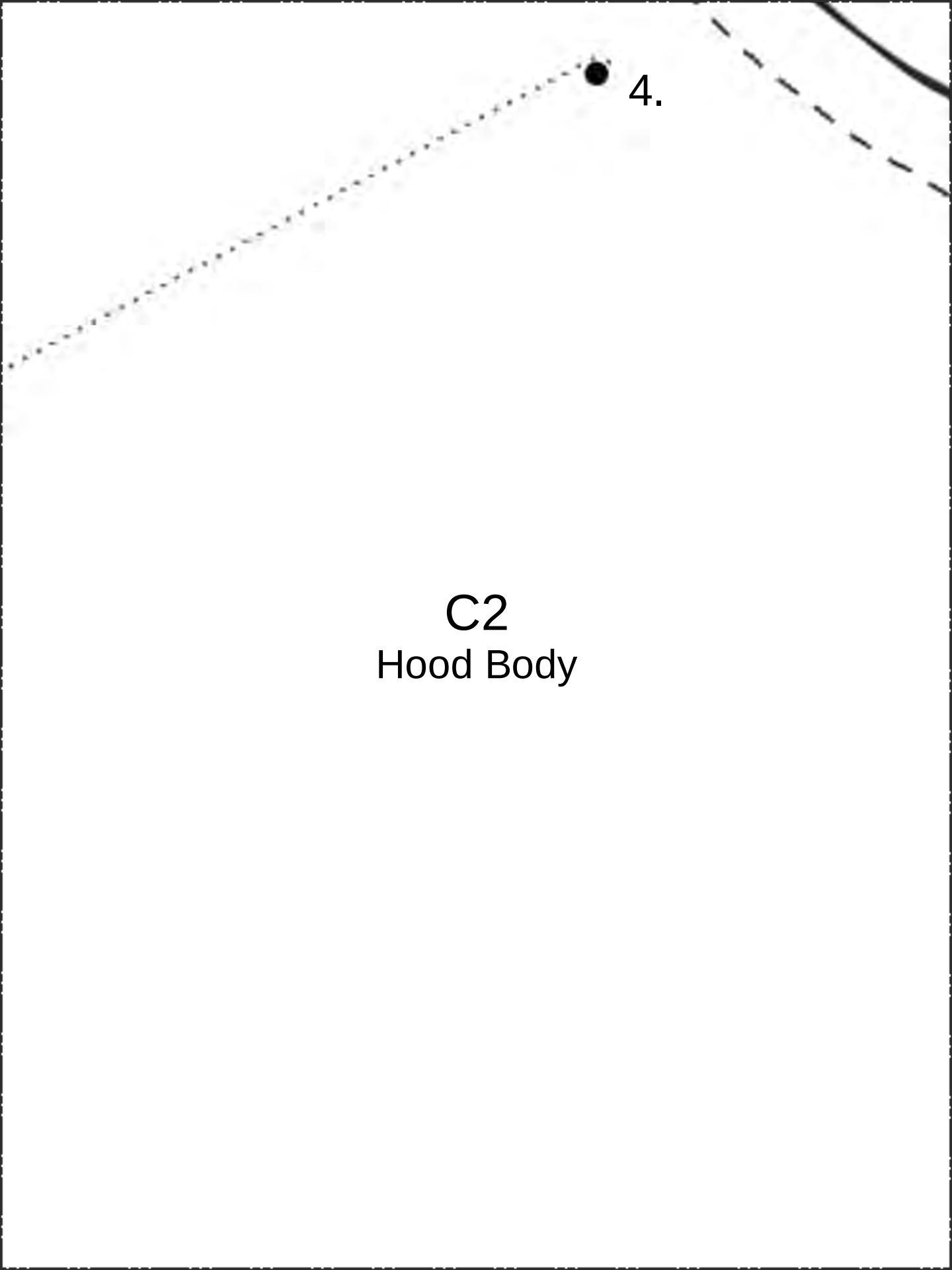


4.


C2
Hood Body

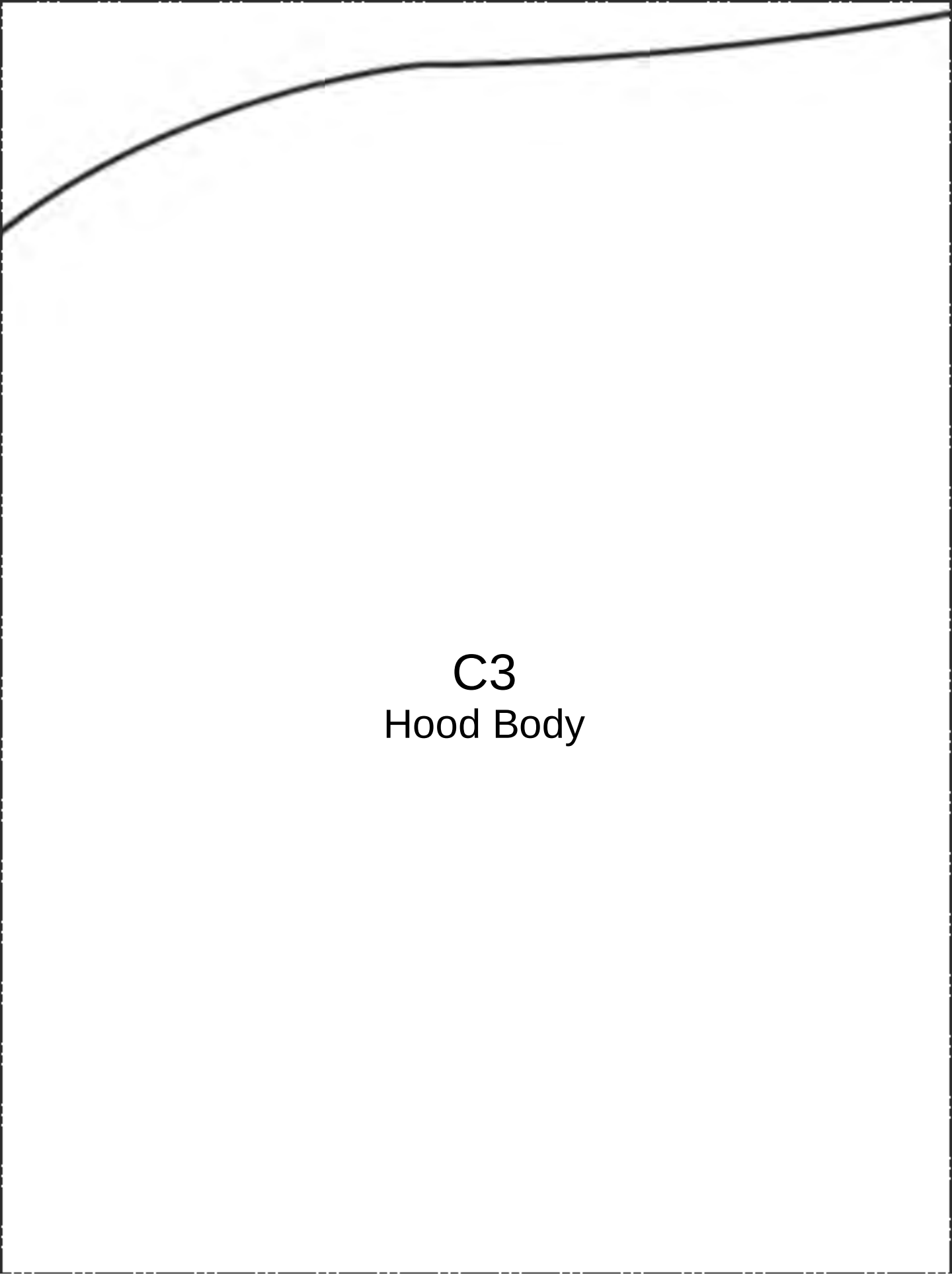

C3
Hood Body


D1 Hood Body

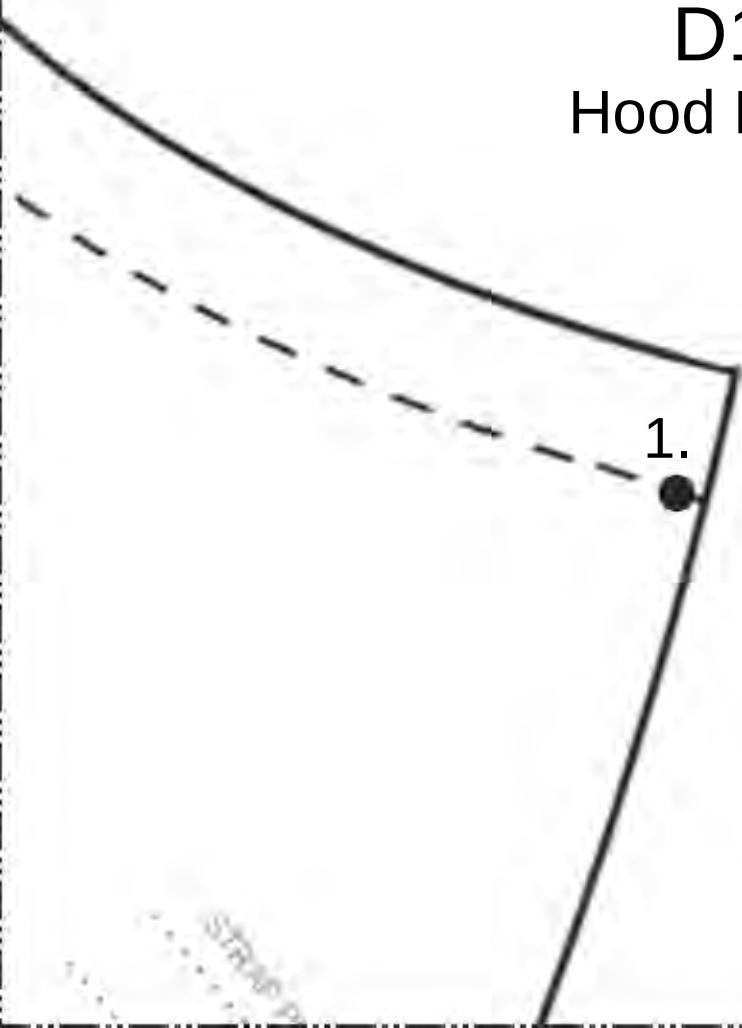

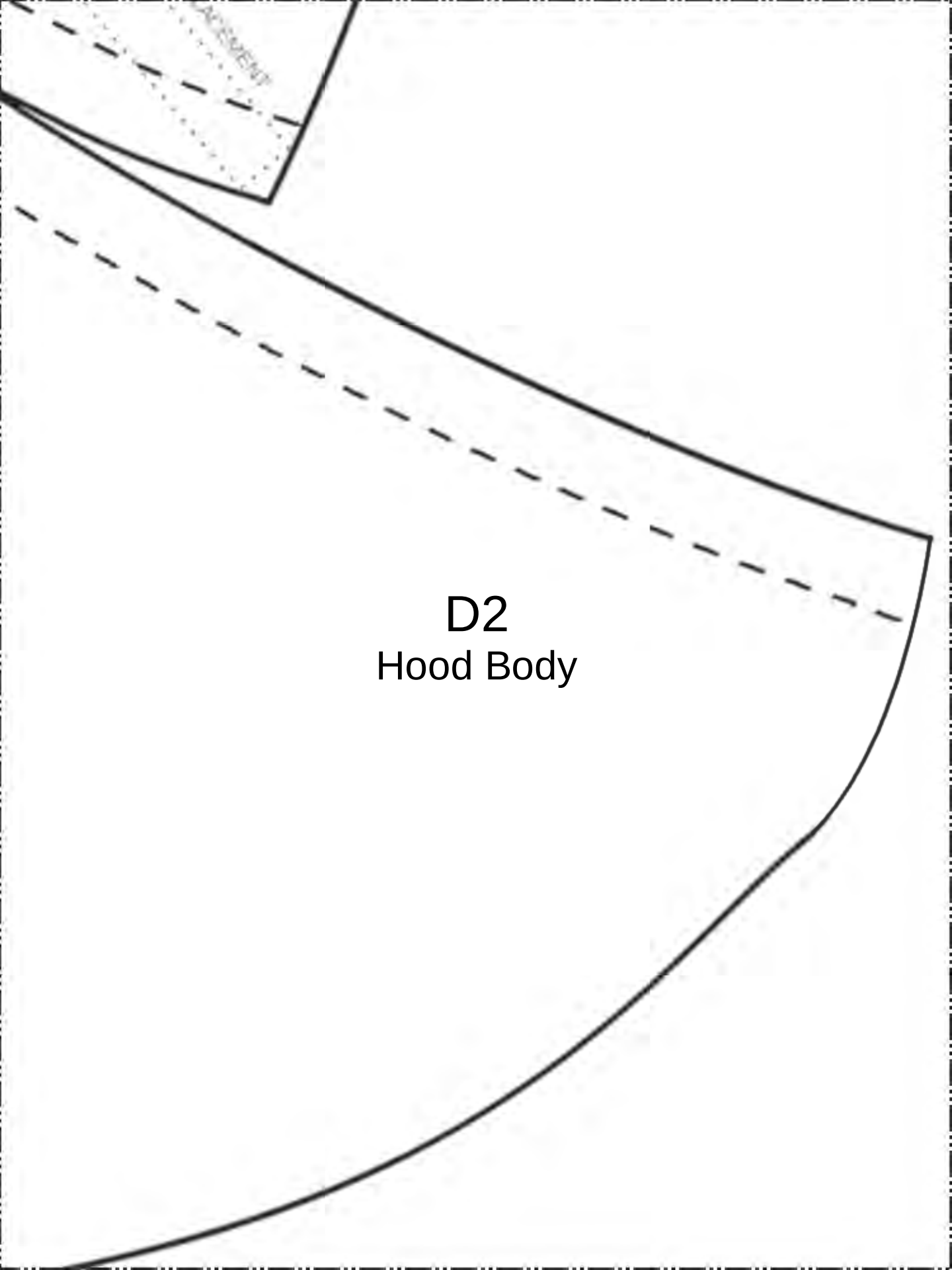


PLACEMENT

D2
Hood Body

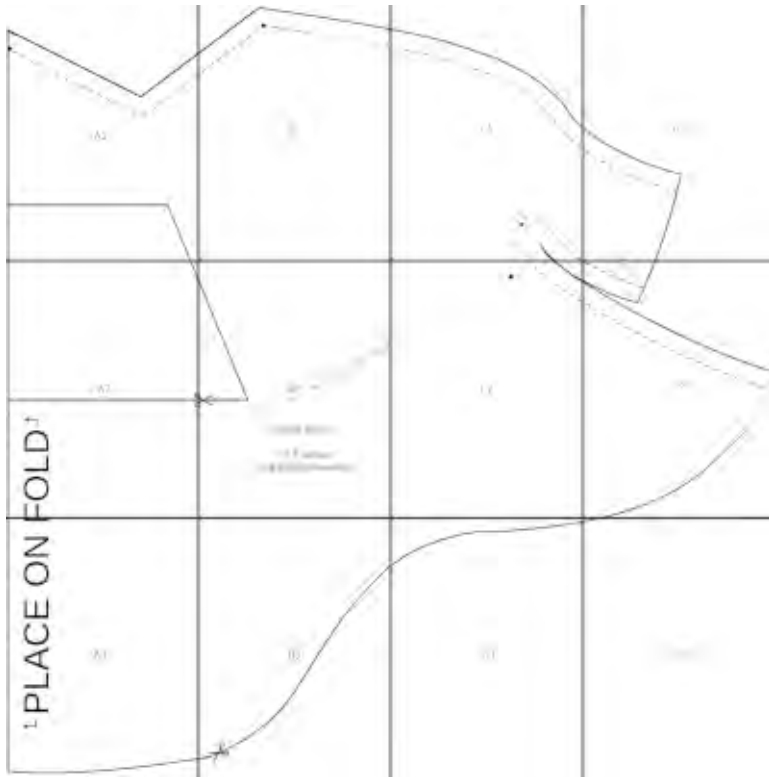

D3

Hood Body

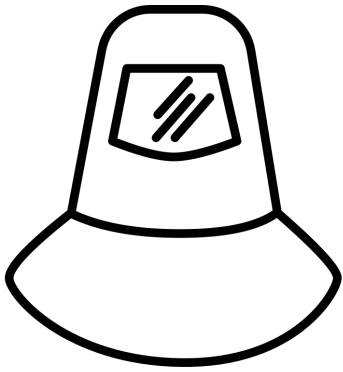

Emergency PAPR Hood

Sanitization Instructions

1. Take out of plastic bag
2. Bake in oven at 70° C (158° F) for 20 minutes

Hoods can be reused as long as they are sanitized before each use

Use Instructions

1. Wear a basecall cap
2. Place hood over head
3. Place forehead strap around head, over top of hat
4. Tuck hood skirt into clothes around neck
5. Insert PAPR hose into trunk
6. Secure hose with rubber band

Created by _____

Optional contact info _____